

# Smart Notes

Thermo Scientific  
CO<sub>2</sub> incubators optimize  
performance by design.



## Which incubation parameters are most important for proper cell growth and expression?

All parameters are important. Proper temperature, gas tension (CO<sub>2</sub>/O<sub>2</sub>/N<sub>2</sub>) and humidity work together to provide optimum growth. Speedy recovery to your set parameters is a critical consideration any time you want to mimic an *in vivo* state.

The goal of placing your cells in a CO<sub>2</sub> incubator at body temperature (37 °C/98.6 °F) is indisputable; mammalian cells will grow best at their native temperature. CO<sub>2</sub> gas serves to maintain *in vivo* pH, similar to CO<sub>2</sub> tension in the bloodstream. High humidity prevents evaporation of growth media. All these parameters work together for healthy cells which express proper protein profiles. Proper culturing is especially important for sensitive primary and stem cells.

But every time you open the incubator door, conditions inside the incubator rush to equilibrate with the conditions in your lab. Once you close the door, the incubator works to re-establish the required state. So how long your cells spend at your specified conditions depends on how long and how often you open the incubator door, and upon the design and engineering of the incubator itself. Different technologies and engineering can result in vastly different recovery times. The Thermo Scientific™ THRIVE™ Active Airflow technology, combined with dual temperature probes, cutting edge CO<sub>2</sub> sensors and a unique covered, integral humidity reservoir, is designed to provide recovery of all parameters within only 10 minutes after a 30 second door opening.



### Different technologies produce different results

Although CO<sub>2</sub> incubators may appear similar from the outside, the engineering inside critically affects your project outcomes. Sensor design and in-chamber location, gas inlet positioning, airflow speed and path, and experienced engineering work together to provide the ideal environment for your cell growth and health.

### Temperature

Cellular functions are highly responsive to temperature. As one example, chick embryos show increased mortality with a difference of only 1 °C<sup>1</sup>. Metabolism and growth slow down at lower temperatures. So fast temperature recovery after a door opening is important, but it's also important to not overshoot the set temperature, because higher temperatures are even more detrimental.

### CO<sub>2</sub> Tension

CO<sub>2</sub> gas works with sodium bicarbonate in the growth medium to control pH to a neutral 7.4. This mimics bloodstream biochemistry. When the pH varies from neutral, cells will first stop growing and then lose viability. Take time in examining your cells, because improper pH may produce morphological changes such as vacuoles in the cytoplasm or granules around the nucleus.

### Humidity

The human body is about 60% water, with internal organs 75-80% water<sup>2</sup>. In a cell culture incubator, balanced growth media provide moisture and nutrients for cells. Humidity of 85-95% limits evaporation of water from the media. Evaporation leaves too-high concentrations of salts, minerals, etc., resulting in toxicity and cell death. High humidity is the most difficult condition to reestablish but is critically important, as evaporation is 4 times faster at 80% humidity than at > 93%<sup>3</sup>.



### Summary

Temperature, gas tension and humidity all work together to replicate *in vivo* conditions. Thermo Scientific THRIVE active airflow, cutting edge in-chamber sensors and design technologies recover all parameters in 10 minutes or less (30 sec. door opening).

[thermoscientific.com/co2](http://thermoscientific.com/co2)

© 2015 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

**Australia** +61 39757 4300  
**Austria** +43 1 801 40 0  
**Belgium** +32 53 73 42 41  
**China** +800 810 5118 or  
 +400 650 5118  
**France** +33 2 2803 2180  
**Germany national toll free** 0800 1 536 376  
**Germany international** +49 6184 90 6000


**India toll free** 1800 22 8374  
**India** +91 22 6716 2200  
**Italy** +39 02 95059 552  
**Japan** +81 3 5826 1616  
**Netherlands** +31 76 579 55 55  
**New Zealand** +64 9 980 6700  
**Nordic/Baltic/CIS countries**  
 +358 9 3291 0200

**Russia** +7 812 703 42 15  
**Spain/Portugal** +34 93 223 09 18  
**Switzerland** +41 44 454 12 22  
**UK/Ireland** +44 870 609 9203  
**USA/Canada** +1 866 984 3766

**Other Asian countries** +852 2885 4613  
**Countries not listed** +49 6184 90 6000

### Figure 1: CO<sub>2</sub> Incubator Design Affects Recovery

Comparison of (A) temperature (within 98% of 37 °C), (B) humidity (to >90%) and (C) CO<sub>2</sub> (within 98% of 5% CO<sub>2</sub>) recovery after a 30 second door opening. Different technologies produce widely varied results. Data from published brochures.

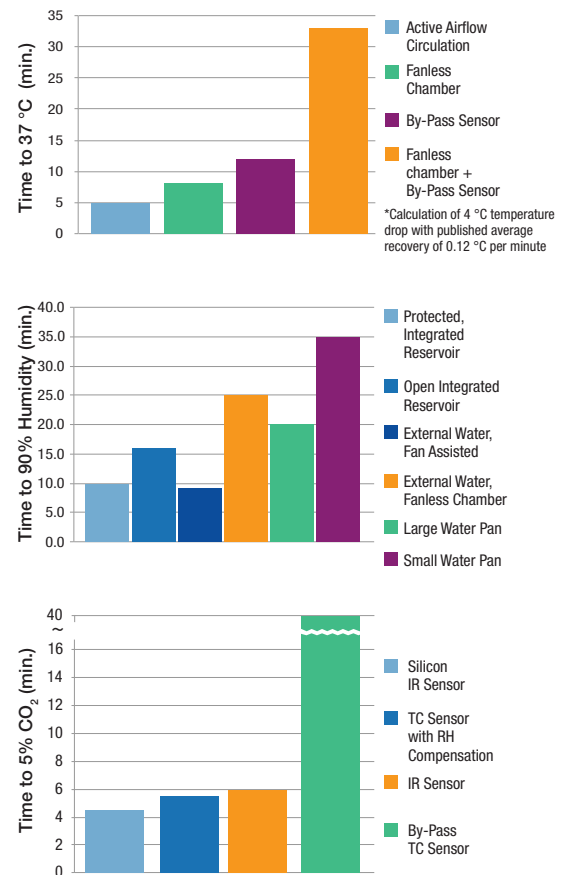

#### References:

1. Poultry Hub: Physiology and Incubation. <http://www.poultryhub.org/most-popular/incubation>
2. Mitchell HH et al. The chemical composition of the adult human body and its bearing on the biochemistry of growth. *J Biol Chem.* 158:625, 1945.
3. Esser P and Weitzmann L. Evaporation from cell culture plates. Thermo Fisher Scientific TILSPNUNCBU02 0111, 2011.